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Research and Production Company
“BIOMEDIS”

**Physiotherapeutic device for application of
bioresonance therapy and anti-parasitic correction of
environment**

“BIOMEDIS M”

TC 9444-001-87530614-2008

Application Instruction

Read instruction before application!

Content

Scope of delivery
Function
Principle of action
Structure of device
Indications for application
Contraindications for application
Side effects
Preparation for usage
Way of application
Recommendations
Safety measures
Technical characteristics
Warranty obligations

Scope of delivery

- 1) Physiotherapeutic device BIOMEDIS
- 2) Network supply battery charger
- 3) USB-cable
- 4) Application instruction (passport)

Functions of device

The medical device “BIOMEDIS M” is indicated for application of bioresonance therapy (BRT), a therapy with the help of impact exercised by electromagnetic oscillations which, under the principle of bioresonance, contribute to restoration of healing and self-regulation’s mechanisms damaged by diseases.

The apparatus “BIOMEDIS M” makes possible to carry out bioresonance exogenous therapy of patient:

- A) With electromagnetic fields, having amplitude modulation of carrying frequency, with high accuracy of setting frequency modulation;

- B) With electromagnetic signals transmitted through proper connector to external interface.

Frequencies of Foll, Schmidt, and Rife are utilized for application of this therapy.

This therapy can be combined with any other, such as pharmacotherapy, physiotherapy, homeopathy, acupuncture, etc.

Each treating program consists of set of frequencies which exercise their impact consistently during specified time duration. Treating complexes are composed of treating programs during programming of the device.

The device may be utilized by several users sequentially.

The device may be equipped with external optional interfaces (such as hand electrodes, plasma generator super-high frequency field).

Principle of device's action and recommendations on application of therapy

Bioresonance therapy represents correction of the body's functions through harmonious transformation of resonance electromagnetic oscillations which suit to living organism for their synchronization.

This treatment is based on suppression of pathologic frequency spectrum oscillations and strengthening of physiological ones, and support of relative synchronization of different wave processes which compose physiologic homeostasis of the body.

The idea of BRT with the help of weak electromagnetic oscillations inherent to patients themselves was expressed for the first time and scientifically explained by F. Morell (1977).

Relative synchronization of different oscillating (wave) processes is kept during normal physiological state, while pathologic states (illnesses) start developing in case of oscillating harmony's deviation in the body.

BRT represents a therapy with electromagnetic oscillations the body's structures come into resonance with. Exposure can be exercised on the level of cell as well as organ, system of organs and the whole body because various levels are characterized and managed by various frequency-wave parameters.

The chief idea of application of resonance in medicine is following: if frequency and form of treating (electromagnetic) exposure are

chosen correctly, normal oscillations can be strengthened in the body, and pathologic ones can be weakened. Thus, bioresonance exposure can be directed for neutralization of pathologic oscillations as well as restoration of physiologic ones damaged during pathologic states; so, hindrances (noises) of informational field of the body can be suppressed: it is impossible to treat any disease without the suppression.

Relative synchronization of various oscillating (wave) processes is supported in body, regardless of age; the body keeps its integrity, balance and health with the help of the processes mentioned above. Imbalance and harmony's deviation of frequency-wave homeostasis appear in case of development of any pathologic pathogenic processes. It causes development of new vicious pathogenic oscillations and resonances which support disease and disturb treating. The program of self-restoration and self-healing given by nature is disturbed.

The task of the device "BIOMEDIS M" is following: to restore correct healthy energetic potential of organs and systems step-by-step, to cause proper resetting, transform unhealthy rhythms and frequencies into healthy row, in accordance with principle of resonance.

The basis of the treating method with the device "BIOMEDIS M" is following: solution of the problem to create a method of electromagnetic therapy which allows causing harmonization in the

work of organs and systems of human body with the help of resonance phenomena caused in the body by specified frequencies.

The essence of the electromagnetic therapy's method is following: restoration of health and liquidation of disease's mechanisms are caused by exposure exercised on the body by impulse low-frequency electromagnetic field with certain carrying frequency modulated by various frequencies.

Success is achieved in treatment because weak electromagnetic field of curative type causes resonance of health and treatment in organs and systems damaged by diseases or because of aging, harmonizing and restoring work of the whole body in general and activity of each cell or organ separately. Algorithm of health replaces algorithm of disease. It is important to notice that it will be your own health: nothing foreign which wouldn't match your body will be implemented from outsides. That's why there is no any chance to overdose or get any side effects because of application of the device. If the frequency sent from the device causes resonance of healing, and if the body needs it, treating effect is caused. If the body doesn't need the frequency, and if resonance doesn't appear, respond doesn't appear too. As a result, no any effect is caused in the body.

The device can be used in house conditions; usage of the device doesn't require contacting a physician or his/her supervision (if

relative contraindications mentioned in guidelines are considered). But if it is needed, contact your personal doctor! Consultation with your doctor can contribute in development of your treating program and increase effectiveness of the therapy.

The therapy is carried out through sessions which compose a course; each session represents a single treating impact of certain program or a complex of programs. Several sessions may be carried out per day. More pronounced effect is caused by a course therapy.

A course therapy includes impact exercised during 14-21 days. A therapy may consist of 5-7 courses with a break of day between therapies, if it is required.

The device should be arranged with its backside to the damaged organ for getting of more pronounced treating effect.

Structure of device

The device represents a generator-radiator of high frequency electromagnetic field of weak intensity with amplitude modulation of carrying frequency. Frequency and form of modulation signal is set through programming (in accordance with chosen algorithm) with the help of a synthesizer. High accuracy of setting modulating frequency is reached with the help of the synthesizer.

The device is equipped with graphical TFT display 1.77 ". The device is manipulated with the help of a keyboard, having tactile feedback.

The device is supplied from built-in LiPol battery with 300 mAh; the battery can be charged through mini USB port of the device, if it is needed. The battery allows the device working 15 hours without stop.

The device has USB port for its connection with PC (if corresponding software presents, USB port is used to download programs and treating complexes). This port can be also used for charging of the device from an external source. A PC or net adapter of 5V, having corresponding connector, can be used as an external source of supply.

The device is equipped with a slot micro SD for extending of its memory. The maximal card memory supported by the device is 8GB.

Indications for application

(in accordance with Methodical recommendations of Ministry Health of Russian Federation №2000/47)

- functional disorders of different genesis;
- diseases of central nervous system and sense organs;
- diseases of vegetative nervous system;

- pain syndromes of various localizations and genesis;
- diseases of cardiovascular system;
- diseases of respiratory system;
- diseases of digestive system;
- diseases of skin and subcutaneous tissues;
- diseases of musculoskeletal system;
- diseases of urinary system and reproduction;
- poorly healing wounds and ulcers;
- other conditions.

Contraindications for application

(in accordance with Methodical recommendations of Ministry Health of Russian Federation №2000/47)

- blood clotting;
- pregnancy (1st trimester);
- benign and cancerous tumors (allowed only under supervision of a highly qualified physician);
- babies that are younger than one year old (administration is allowed only under supervision of a physician);
- existence of transplanted organs (allowed under supervision of a doctor);
- existence of the implanted electrocardiostimulator;
- individual intolerance to electric current (if hand electrodes are used);
- epilepsy (episode), convulsive attack;
- congenital defects of central nervous system;
- severe psychic disorder or intoxication.

Reactions of body to process of therapy

(mainly for anti-parasitic correction)

Aggravation of the main disease can take place during the therapy; the exacerbation may be accompanied with general malaise, weakness, diarrhea, frequent urination, thirst, sleepiness, temperature reaction, skin rashes, etc. It is needed to prolong interval between sessions for improvement of physical condition in this case.

The all reactions may be observed only during the first two weeks of application of the device "BIOMEDIS M". There is no any necessity to use additional medications under this circumstance.

These reactions are connected with individual response of the body to wave impact of the device, and manifest increase of excretory reactions of the polluted body.

Sometimes reactions may appear later and discomfort condition may be presenting during a month after you have begun to use the device.

ATTENTION!

Consider contraindications listed in the guide for application of the device! Contact your doctor!

Preparation for usage

Image 1. Front panel of device

- 1 - On/off button
- 2 – Display
- 3 – Arrow keys

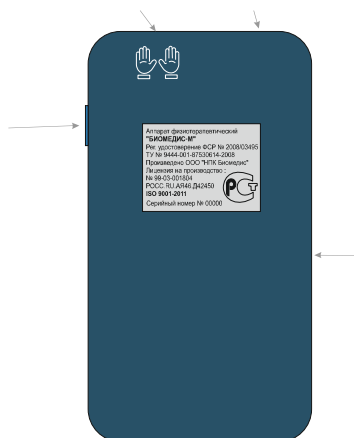
Image 2. Rear panel of device

- 1 – On/off button
- 2 – Connector of external interface
- 3 – mini USB connector
- 4 – microSD connector

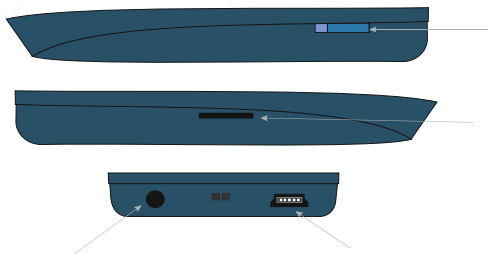
Image 3. Side panels

- 1 – On/off button
- 2 – microSD connector
- 3 – Connector of external interface
- 4 – mini USB connector

2. Check workability of the device, turning on power! Move the button (Images 1, and 2) to the ON position for this purpose! When you do it, the name of the device will be written on display. Then you will see the MAIN MENU (Image 4) of program on the display of the device. If it is needed, charge the device (see the paragraph 6)!

Move the button 1 to the position OFF (Images 1, and 2) in order to turn off the device!

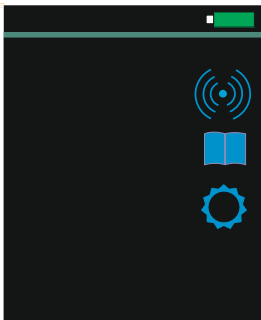
Image 4. Main Menu of device

Way of application

12. Main Menu

12.15. After you will have turned on the device, the following information will be displayed on the screen (from the top till the bottom) (Image 4):

Therapy - Regime of Therapy

T-Book - Regime of viewing text description of program

Setting – Regime of setting device

12.16. Use  and , and then click

for choice of the needed item!


Image 5. Choice of memory's directory

13. Regime Therapy

13.15. If a memory card is inserted into the device, after you will have chosen the regime Therapy, you will be asked which memory's directory should be used for reading files: the directory of the device itself – Local Complex list, or the directory of micro SD – Card Complex list (Image 5). Choose needed but-

ton  or , and then press ! If one of the directories is empty, this item will be not contained in the menu. If you desire to return to the menu, you should choose “back to menu” and then click .

13.16. Then Complex list will be displayed in the chosen directory (Image 6).

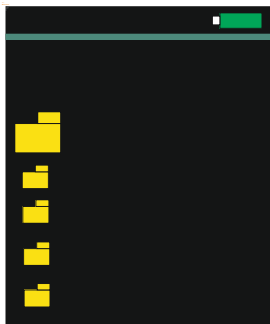
Image 6. Display of Complex list

13.17. Use the buttons  and , and then press  for choice of the needed complex!

13.18. Click “back to menu”, and then press  to return to menu!

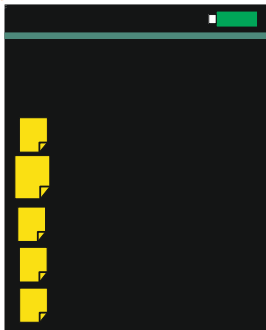
Image 6.1 Display of complex contain (list of programs contained in the chosen complex)

13.19. Choose a program which will be used as a start of the

complex, utilizing the buttons

and

and then press the button

(Image 6.1)!

13.20. If you want to go back to the previous menu item, choose the top menu item, having the name of the complex, and then

press the button !

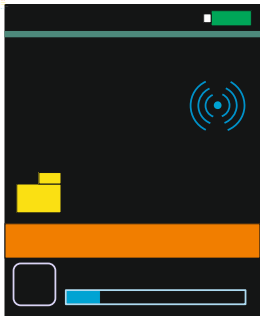
Image 7. Display of complex's time parameters

13.21. After you will have chosen a program to start a complex, the following information will be displayed on the screen of the device from the very top till the bottom (Image 7):

- time duration left till the end of the complex in the format "hh:mm:ss";
- time duration left till the end of program in the format: "hh:mm:ss";

- name of the complex applied;
- name of the program applied;
- regime of the work (pause-application);
- time duration of the program applied with graphical display of fine-tuning. The program is not

carried out in this regime. Click  to go back to the previous menu item!

13.22. Press the button  (picture 8) to launch the complex (program)!

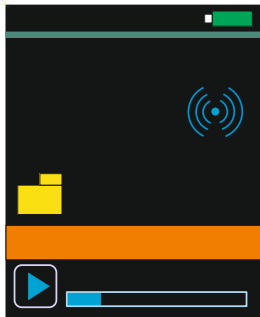
Image 8. Implementation of complex

13.23. The pictogram  in the left bottom of the display is

transformed to the pictogram  during implementation of complex. Timer counts down duration of complex (program).

13.24. Press  to pause implementation of complex (pro-

gram)! Press , and then click  to stop complex (program) completely!

13.25. Use the buttons   (Image 9) for fast transition to previous/next program during implementation of complex!

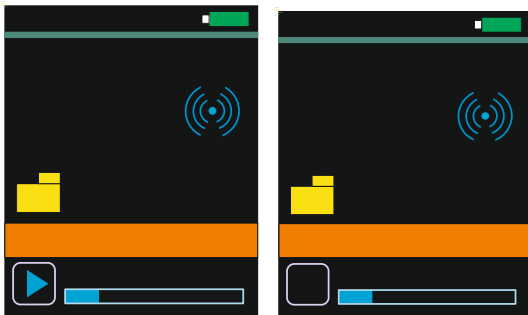
Image 9. Changing of programs inside complex

13.26. Complex (program) can be implemented through 3 regimes; each regime is displayed through corresponding pictogram on the top of the display next to the record "Therapy":

- normal regime; complex is implemented once consistently: one program after another;

- this pictogram shows that just one program included in complex, displayed on the screen, is implemented cyclically till the

forced stopping. Choose needed program, using

the buttons  and !

-  - cyclical regime: complex is implemented multiply: one program after another, till the forced stopping.

Press  during implementation of complex (program) to choose the needed regime for implementation of complex (program) (Image 10)!


Image 10. Choice of regime for implementation of complex

Choose the upper menu item, and press , to go back!

2.13 Choose "Repeat" and press  (Image 11), to continue setting!


Image 11. Choice of regime for implementation of complex

2.15 Then choose needed regime:

- Normal – normal regime
- Repeat one – repeat of one program

- Complex – repeat of the whole complex

through usage of the buttons

and pressing of

- 2.16 Choose the upper menu item and press

to go back to implementation of complex (program)!

14. Regime T-Book

- 14.15. This regime allows viewing files – descriptions of each therapeutic program. Files are created automatically with the help of a special program installed on PC. Choose T-Book and

press  (Image 12)!

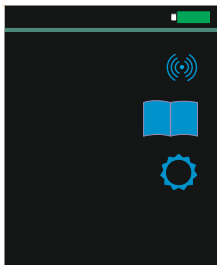
Image 12. Choice of regime T-Book

If it is needed, repeat instructions given in the paragraph 2.1!

14.16. Then the complex list is displayed in the chosen directory (Image 13).

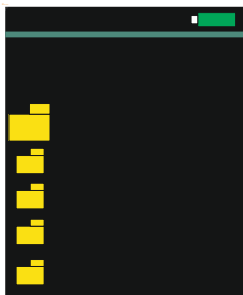
Image 13. Display of complex list

14.17. Use the buttons   and press

 for choice of the needed complex!

14.18. Choose “back to menu” and press  to return to menu!

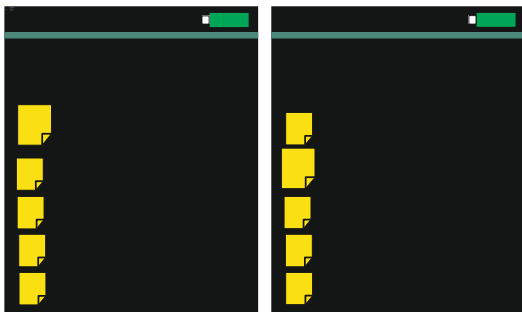
Image 14. Display of complex contain (list of programs included in chosen complex)

14.19. If you desire to view description of a program, choose the

program, using the buttons  , and then

press  (Image 14)!

14.20. If you want to return to previous menu item, choose the upper menu item, having the name of the complex, and then

press !


Image 15. Display of text description of program

14.21. Then text description of program will be displayed (Image

15). Use the buttons   to scroll through the text!

15. Regime Setting. Settings

15.15. This regime gives the possibility to change settings of the device. Choose the item “Setting” in main menu, and then

press  (Image 16)!

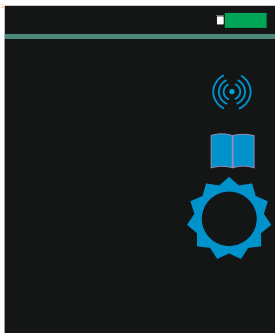
Image 16. Choice of regime Setting

15.16. The list of settings includes the following items (Image 17):


Image 17. List of settings

«LCD Set» – setting of time, after which keyboard will be blocked and display will be switched off

«Language» – menu language

«Memory info» – information about memory

«Firmware version» – firmware version of device

«Firmware upgrade» – updating of firmware (only for service centers).

Choose needed item, using   and pressing

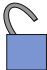
! Choose “back to menu” and press  to go back to main menu!

- 15.17. **LCD Set** is time duration, after which the keyboard will be blocked and the display will be switched off (for charged energy saving). Possible time options last from 0 to 50 seconds. If you set the time 0, the display and the keyboard are not

blocked. Use the buttons  , and press

 to set up needed option! After the display has been switched off and the keyboard has been blocked during working of the device, pressing of any button causes activation of the display; then if you try to use the keyboard, the sign of

locked padlock  appears on the display. Press and hold

 till the sign of opened padlock  appears, in order to unblock the display.

- 15.18. **Language** is menu language. Use the buttons

 , and press  to set up needed language!

15.19. **Memory info** displays memory capacity of the device, noticing used and unused memory.

15.20. **Firmware version** displays the firmware's version of the device.

16. Choice of impact's intensity

16.15. This version of the device allows the possibility to change intensity of the device's high-frequency radiation. Intensity of the device can be set up, beginning from 0 till 10 nominal units. The level 10 is used as default in the device. To change the level of radiation during the process of the therapy's appli-

cation, you should press , select the needed level,

using the buttons  , and then press

.

17. Charging of device

17.15. The pictogram which displays the level of battery's charging is in the upper right corner. If the battery's intensity drops below the minimum permissible level, the pictogram is EMPTY. It is needed to charge the battery in order to avoid switching off the device during the therapy.

Attention! Starting charging can be carried out without interrupting the process of the therapy, only if you use battery charger. In this case you should use the charger, utilizing the corresponding device's connector with the help of the cable included in the kit!

17.16. The device has USB-port for connection with PC. The port can be used as an input for charging of the device from external supply source and charging of the device's battery. PC or any adapter of 5V, having appropriate connector, can serve as external supply source. If the device is connected with PC, the record CONNECT is displayed on the screen.

(Attention! If the therapy is being carried out, the process of the therapy will be stopped, you will need to restart the process of the therapy after prolonged holding and releasing

the button  and going to main menu.)

- 17.17. PC determines the device as standard external drive – drives (depending on the condition, whether SD-card is inserted into the device or not).
- 17.18. The process of charging is displayed through pulsation (filling) of the battery icon. When the battery is charged, the icon of battery capacity will be full.

Recommendations

1. Try to avoid serious physical activities on the days of sessions!
2. Regardless of time duration taken by sessions per day, drink not less than 1 liter of clean not boiled water a day! It is recommended consuming 30ml of water per 1kg of body mass.
3. If you have digestive problems (diarrhea or constipations), you should take corresponding biologically active additives, tea remedies and products which contain lacto- bifido-bacteria under recommendations of your treating doctor.
4. If you have acute diseases (poisoning, ARI, flu), it is recommended applying sessions in each 2 hours till you get improvement and achieve relief of diseases' symptoms.
5. Application of the device "BIOMEDIS M" doesn't replace usage of medical therapy and remedies of natural medicine; on the contrary, it increases their effectiveness.

6. Since intensity on electrodes is quite low, application of sessions doesn't cause any feelings, or it causes just slight tingling. Arrange cosmetic cotton pads soaked in water under electrodes for increase of treatment's effectiveness!
7. It is recommended to use cosmetic cotton pads soaked in water arranged on locations of electrodes' contacts, if skin damages present.
8. However hands of treated person are not engaged during the therapy, it is recommended to rest and turned himself/herself to the process of healing.
9. Needed quantity of sessions is chosen, basing upon individual tolerance.
10. The device may be used after vaccination, only after contacting physician.
11. In case of diseases which are dangerous for health, the device may be applied only under supervision of health care provider.

Conclusion

Success can be achieved by usage of the device "BIOMEDIS M" in case of any disease which doesn't require urgent surgical intervention, if person has been diagnosed correctly. It is necessary to contact doctor to use the device more effectively.

Indications for safety measures

12. The device has been made in accordance with requirements for devices, having internal source of supply, according to State Standard P502 67.0-92. "Way of protection from damage by electric current".
13. Water must not penetrate into inputs and connectors of the device.

14. In case of improper work of the device, you should turn off it immediately. Don't use devices with have evident damages (broken wire or broken casing)!
15. Don't keep the device near open flame!
16. Don't leave the device prolonged time under direct sun lights!
17. Keep the device away from children!
18. Don't use the device in bathroom or in another wet room!
19. The device can be used for children's treatment only under supervision of adults.
20. Don't put the device into ovens of SHF!
21. Don't disassemble the device!
22. Don't put the device into water!

Technical characteristics

Dimensions (not over than) 79x42x11 mm

Parameters of therapy's signal on connector of external interface:

Amplitude of signal (not over than) 3V

Form of signal - any

Frequency range - 0.01 - 20000.00 Hz with frequency step 0.01 Hz

Parameters of signal for contactless therapy:

Carrying frequency - 24 MHz+- 1%;

Intensity of electromagnetic field while being not over than 10 cm from the device - 10 V/m;

Type of modulation - amplitude modulation;

Frequency range of modulating signals 0.01-20000.00 Hz with frequency step 0.01 Hz;

Form of modulating signals - any;

Supply voltage 3.7 V (inner LiPol battery with capacity 300 mAh);

Current consumption - not over than 30 mA;

Memory capacity of device - 8 Gb;

Maximal memory of microSD card supported by device - 8 Gb.

All needed connecting cables come in the kit together with the device.

Warranty obligations

6. The manufacturer guarantees conformity of the device with technical conditions, if conditions of storage, transportation and exploitation are followed.
7. Warranty period of operation is one year since the day of purchase.
8. If mechanical failures of components of the device "BIOMEDIS M" occur through no fault of user, replacement of the whole device as well as its components is possible.
9. The warranty doesn't apply to device which has signs of essential mechanical and chemical damages which destroy constructive integrity of the device.
10. Warranty period of operation is calculated since the date of purchasing device or since the date of getting the device by postal service.

Acceptance report

The device "**BIOMEDIS M**" has been manufactured and accepted in accordance with regulatory requirements of state standards and acting documentation; and it is considered to be operable.

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